

Variolink Esthetic® LC


Placement of inlays (thickness < 2 mm) (glass-ceramic, e.g. IPS e.max®)

Try-in



Variolink
Esthetic
Try-In





Restoration

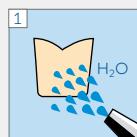



Monobond
Etch & Prime

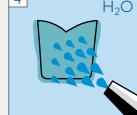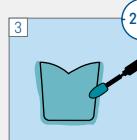

Preparation

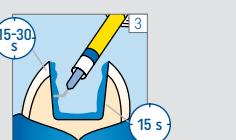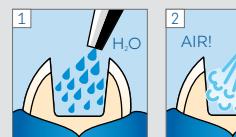

Total Etch:
optional

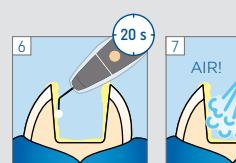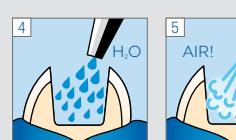


Adhere
Universal

Avoid pooling


Bluephase
G4

Cementation

Variolink
Esthetic LC

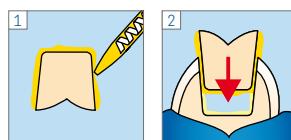


Bluephase
G4

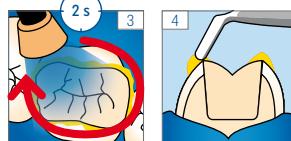


Secure the restoration until final polymerisation

Liquid
Strip

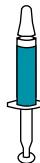

* per segment and mm ceramic

Fluor
Protector S

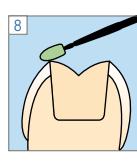

Manufacturer
Ivoclar Vivadent AG
Bendererstrasse 2
9494 Schaan/Liechtenstein
www.ivoclar.com



See Instructions

Date information prepared:
2022-08-23
681584/EN/Rev. 1



more information
www.cementation-navigation.com

ivoclar

